

# DC Series

## Fluid Warming Cabinet

### Temperature range

32°-40°C (90°-104°F)  
adjustable injection mode

32°-66°C (90°-150°F)  
adjustable irrigation mode

### WarmRight® technology

uses multiple zone heating  
(Patent number US 8,217,316)  
that monitors surfaces for  
safe temperature ranges

Efficient, balanced heating  
system allows for low  
energy consumption  
and minimal heat loss

Interior lighting casts  
a comforting blue glow  
with two (2) different  
intensity settings

Ultra quiet air flow  
system designed for  
balanced air movement

Available in two (2)  
different cabinet sizes

Available with two (2)  
baskets (standard),  
one (1) basket (optional), or  
without basket (optional)



**DC400L**

**DC250L**

(SHOWN WITH OPTIONAL  
ELECTRONIC COMBINATION DOOR  
LOCK AND SINGLE BASKET INSERT)

 **ENTHERMICS**  
**Medical Systems**

An ISO 13485:2003 certified company

*the warming people®*  
FLUID & BLANKET WARMERS

[www.enthermics.com](http://www.enthermics.com)

*Experience the Entermics difference:  
Quality beyond comfort.*

# DC Series Fluid Warming Cabinet

## Specifications:

- Warming cabinet with stainless steel interior
- Ultra quiet air flow system designed for balanced air movement
- Energy efficient, e-coated glass window in door allows for inventory observation
- Easy, hands-free, push-button door release
- Door is fully gasketed and hinged on the right side of the unit
- WarmRight® technology employs a multi-zone heating system that monitors cavity surfaces to ensure a safe temperature range
- Furnished with four (4) 3/4" (19mm) non-skid feet

## Control:

- The warming cabinet can be set to warm either irrigation fluids (IRR) or injection fluids (INJ), with separate temperature ranges depending on the choice selected. IRR temperature may be adjusted from 32° to 66°C (90° to 150°F), and the INJ temperature can be adjusted from 32° to 40°C (90° to 104°F)
- Efficient, balanced heating system allows for low energy consumption and minimal heat loss
- Operates in Celsius or Fahrenheit
- Four digit L.E.D. display
- Actual temperature indication
- Built-in speaker for audible feedback
- User Interface lock-out control feature

## Additional features:

- L.E.D. interior lighting casts a comforting blue glow with two (2) different intensity settings, and off mode.
- Convenience: In the event of a power failure the cabinet will remember its settings and begin to operate as before when power is restored with a visible indicator.
- Safety: A thermal shut-off system, separate from the electronic control, is included as an additional safety feature.
- Stackable design - see options and accessories list for available combinations.

## Factory Installed Options

### DC250L

- ☐ **Casters, 3" (76mm)** [5012693]  
OVERALL HEIGHT INCREASES 3.6" (91.4mm)

**Door Lock** – ALLOWS THE USER TO LOCK THE DOOR.  
*AVAILABLE WITH ORIGINAL ORDER ONLY.*


- ☐ **Electronic combination lock**  
(INCLUDES MANUAL OVERRIDE) [5015173]
- ☐ **Manual key lock** [5016065]
- ☐ **Solid Door Assembly** [5015697]
- Stacking Kits** (OTHER OPTIONS AVAILABLE, CONTACT FACTORY)
  - ☐ **DC250L over DC250L or DC350** [5013021]
- ☐ **Stainless Steel Exterior**

### DC400L

- ☐ **Casters, 3" (76mm)** [5012693]  
OVERALL HEIGHT INCREASES 3.6" (91.4mm)

**Door Lock** – ALLOWS THE USER TO LOCK THE DOOR.  
*AVAILABLE WITH ORIGINAL ORDER ONLY.*


- ☐ **Electronic combination lock**  
(INCLUDES MANUAL OVERRIDE) [5015173]
- ☐ **Manual key lock** [5016065]
- ☐ **Solid Door Assembly** [5015698]
- Stacking Kits** (OTHER OPTIONS AVAILABLE, CONTACT FACTORY)
  - ☐ **DC400L over DC400L or DC750** [5013024]
- ☐ **Stainless Steel Exterior**

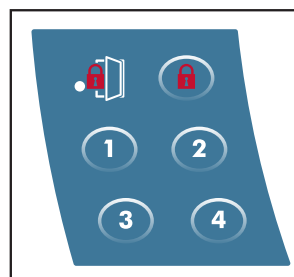

Electronic Combination  
Door Lock


# DC Series Fluid Warming Cabinet

## DC250L

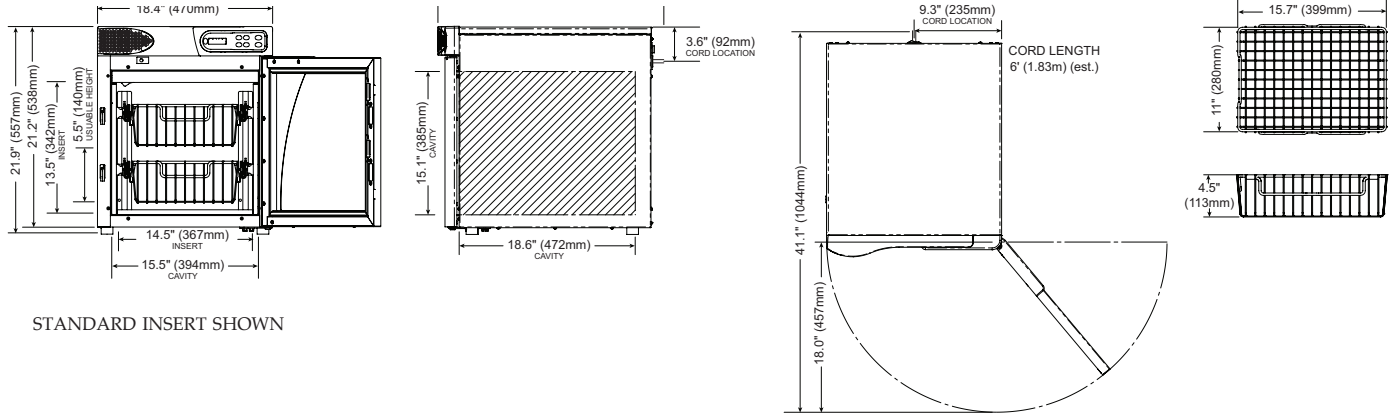

STANDARD INSERT SHOWN

## DC400L

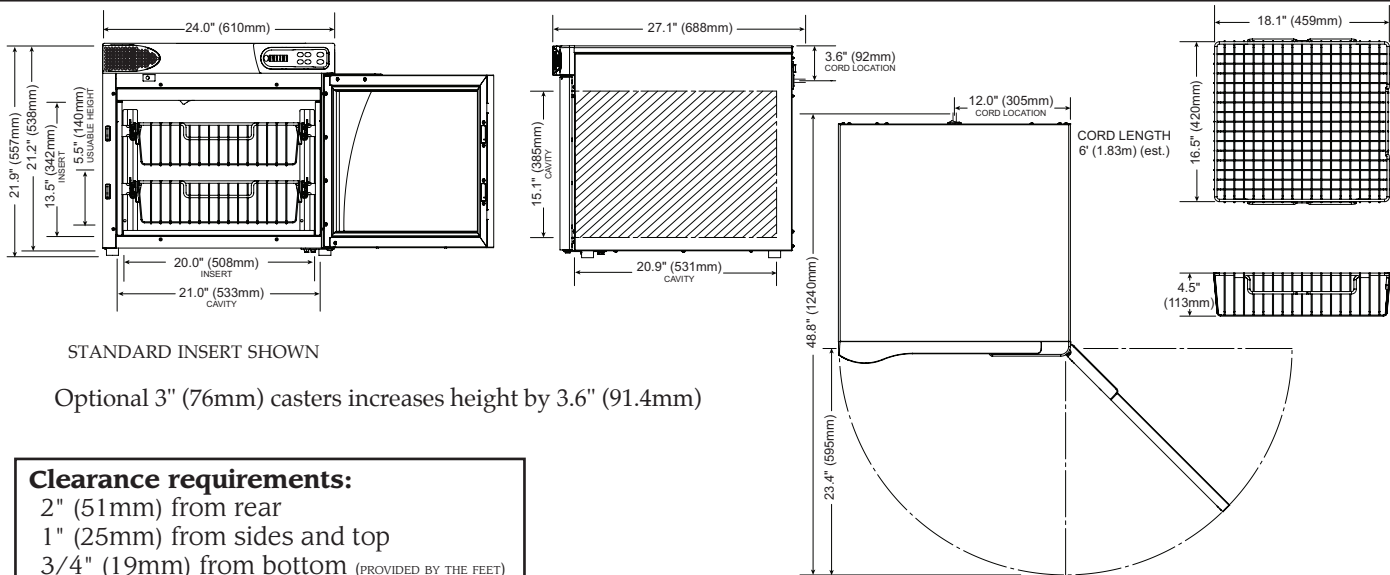

STANDARD INSERT SHOWN

Optional 3" (76mm) casters increases height by 3.6" (91.4mm)

### Clearance requirements:

2" (51mm) from rear  
1" (25mm) from sides and top  
3/4" (19mm) from bottom (PROVIDED BY THE FEET)

### DC250L

### DC400L

Dimensions (H x W x D)	21.9" x 18.4" x 24.8" (557mm x 470mm x 630mm)		21.9" x 24.0" x 27.1" (557mm x 610mm x 688mm)	
Fluid Capacity	Two (2) Baskets, Standard	8 (1) liter bags per basket/16 bags total	Two (2) Baskets, Standard	14 (1) liter bags per basket/28 bags total
	One (1) Basket (POSITIONED ON BOTTOM RAIL), Option	12 (1) liter bottles	One (1) Basket (POSITIONED ON BOTTOM RAIL), Option	20 (1) liter bottles
	Without Basket, Option	CONTACT FACTORY	Without Basket, Option	CONTACT FACTORY
Weight	Net: 85 lb (39 kg)	Ship**: 135 lb (61 kg)	Net: 100 lb (45 kg)	Ship**: 165 lb (75 kg)
Electrical	120 V.A.C. — 60 Hz, 1 ph 0.7 kW, 5.8 Amps Safety Class I Equipment No Applied Parts Mode of Operation: Continuous	NEMA 5-15P 15A - 125V Plug Hospital Grade IP X0	120 V.A.C. — 60 Hz, 1 ph 0.7 kW, 5.8 Amps Safety Class I Equipment No Applied Parts Mode of Operation: Continuous	NEMA 5-15P 15A - 125V Plug Hospital Grade IP X0
	220 V.A.C. — 50 Hz, 1 ph 0.6 kW, 2.7 Amps Type B Equipment	CEE 7/7* 220-230V Plug	220 V.A.C. — 50 Hz, 1 ph 0.6 kW, 2.7 Amps Type B Equipment	CEE 7/7* 220-230V Plug
	230 V.A.C. — 50 Hz, 1 ph 0.7 kW, 3.0 Amps Type B Equipment	BS 1363 Plug* (UK only) CEE 7/7* 220-230V Plug	230 V.A.C. — 50 Hz, 1 ph 0.7 kW, 3.0 Amps Type B Equipment	BS 1363 Plug* (UK only) CEE 7/7* 220-230V Plug

\*Other international plugs are available, contact factory for more information.

\*\*Domestic ground shipping information. Contact factory for export weight and dimensions.



**ENTHERMICS  
Medical Systems**

An ISO 13485:2003 certified company

**Enthermics Medical Systems - The Warming Company®**  
PO Box 443 · Menomonee Falls WI 53052-0443  
W164 N9221 Water St · Menomonee Falls WI 53051  
Phone 262-251-8356 · 800-862-9276 · Fax 262-251-7067  
generalinfo@enthermics.com  
www.enthermics.com · www.ivnow.com