

ivNow®

Patent WO 2010/030746
Canada Patent 2692640

Modular configuration

Wall mount: 1-6 bag capacity

Countertop: 1-3 bag capacity

Pole mount: 1-3 bag capacity

Accommodates 0.5L, 1L, 2L & 3L bags

Anodized aluminum heat plate

Lightweight, medical grade plastic housing

Intelligent bag sensors

Dual temperature sensors

Heats to 40°C (104°F) in 30 minutes or less*

*3L bags require slightly longer heating time

Accuracy of +0°/-2°C (+0°/-3°F)



 **ENTHERMICS**
Medical Systems


www.ivnow.com

*Experience the Enthermics difference:
Quality beyond comfort.*

ivNow® Modular Fluid Warmer

The ivNow fluid warmer quickly warms and maintains the temperature of injection/intravenous and irrigation solutions prior to their use. The specially contoured warming module cradles solution bags in 0.5, 1, 2 & 3 liter sizes. 3-liter bag capacity available on ivNow-1, ivNow-2, & ivNow-3 with an additional adapter. Individual ivNow units cannot be stacked in the field. Multiple cavity units are available from the factory in the configurations below. The unit is controlled by one (1) power switch and individual electronic controls with L.E.D. display for each cavity. The control can easily be set to display temperatures in Celsius or Fahrenheit. A sensor in the heating

plate detects the presence of a bag and engages the heating mechanism to quickly begin warming the fluid. Two (2) temperature sensors work in unison to precisely and continuously read the temperature of the bag and another sensor monitors the plate temperature. A green ready light will illuminate when fluid is within $\pm 0/-2^{\circ}\text{C}$ ($\pm 0/-3^{\circ}\text{F}$) of the set point temperature. The heater will reengage as necessary to maintain the temperature within $\pm 0/-2^{\circ}\text{C}$ ($\pm 0/-3^{\circ}\text{F}$) of the set point. The electronic control monitors the length of time the bag has been held at temperature, beginning when the bag reaches set point temperature. A status button will display the time the fluid has been held at temperature.

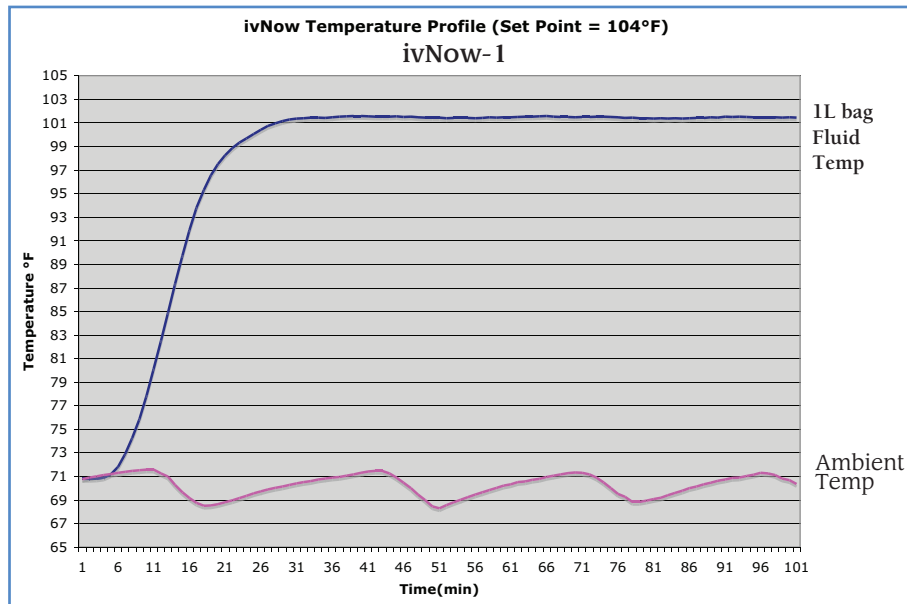

Available Configurations

Note: Individual ivNow units cannot be stacked in the field. Multiple cavity units are available from the factory in the configurations below. 3-liter bag capacity available on ivNow-1, ivNow-2, & ivNow-3 with an additional adapter.

ivNow-1 (one (1) bag capacity) can be placed directly on a countertop, mounted on an optional countertop pedestal stand, mounted on a mobile equipment pole stand, or mounted on a wall using brackets*.

ivNow-2 (two (2) bag capacity) can be placed directly on a countertop, mounted on a countertop pedestal stand, or mounted on a wall using mounting brackets*.

ivNow-3 (three (3) bag capacity) must be mounted on a countertop pedestal stand, a heavy duty mobile equipment pole stand, or on a wall using mounting brackets*.

ivNow-4 (four (4) bag capacity) must be mounted on a wall using mounting brackets*.

ivNow-5 (five (5) bag capacity) must be mounted on a wall using mounting brackets*.

ivNow-6 (six (6) bag capacity) must be mounted on a wall using mounting brackets*.

*Mounting hardware sold separately.


ivNow® Modular Fluid Warmer

Specs	ivNow-1	ivNow-2	ivNow-3	ivNow-4	ivNow-5	ivNow-6	Plug	Certification
Height (wall unit)	8.0" (203.1mm)	13.4" (340.5mm)	18.1" (477.9mm)	24.2" (615.3mm)	29.6" (752.7mm)	35" (890.1mm)	—	—
Width (wall unit)	14.5" (368.7mm)	14.5" (368.7mm)	14.5" (368.7mm)	14.5" (368.7mm)	14.5" (368.7mm)	14.5" (368.7mm)	—	—
Depth (wall unit)	6.74" (171.2mm)	6.74" (171.2mm)	6.74" (171.2mm)	6.74" (171.2mm)	6.74" (171.2mm)	6.74" (171.2mm)	—	—
Cavities	1	2	3	4	5	6	—	—
Net Weight	8 lbs (4 kg)	14 lbs (6 kg)	20 lbs (9 kg)	25 lbs (11 kg)	31 lbs (14 kg)	37 lbs (17 kg)	—	—
Ship Weight	11 lbs (5 kg)	16 lbs (7 kg)	23 lbs (10 kg)	28 lbs (13 kg)	34 lbs (16 kg)	40 lbs (18 kg)	—	—
120V Amps	1.3	2.5	3.8	5.0	6.3	7.5	NEMA 5-15P 15A - 125V Plug Hospital Grade	UL CSA US
120V Watts	150	300	450	600	750	900		
230V Amps	1.7	1.3	2.0	2.6	3.2	3.9	BS 1363 Plug (UK only) CEE 7/7 220-230V Plug	CE
230V Watts	150	300	450	600	750	900		

Note:

- Countertop unit (1 or 2 cavity), power cord exits back of unit.
- Wall mounted units (up to 6 cavities), power cord exits bottom of unit (shown below).
- Mobile equipment pole mounted unit (up to 3 cavities), power cord exits bottom.
- Countertop pedestal mounted units (up to 3 cavities), power cord exits back.

* Cord exits back, maintain a minimum of 4" (102mm) clearance for bending of cord.

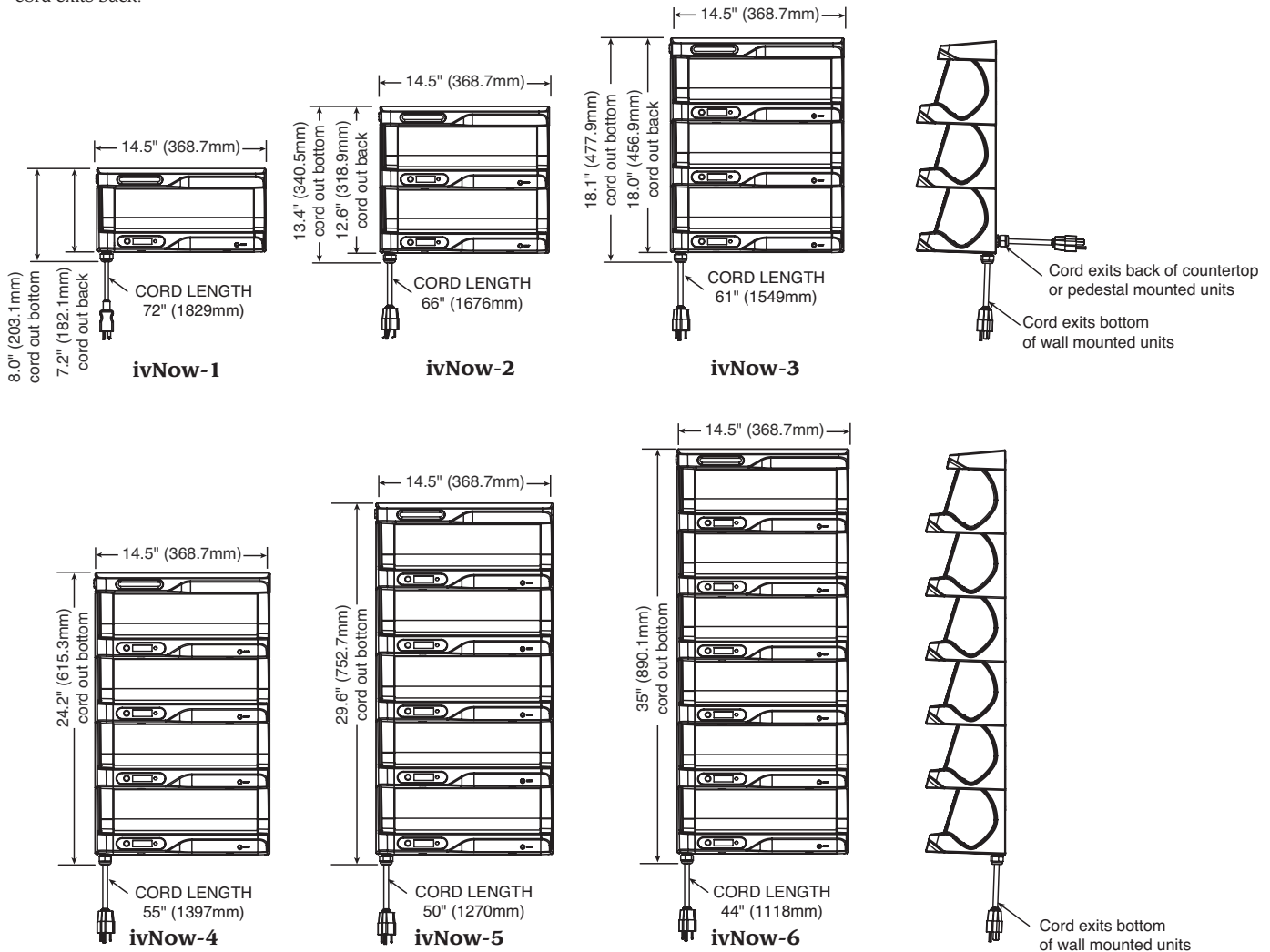

NOVATECH
MEDICAL

1380 Joliot-Curie, bureau 808
Boucherville (Québec) J4B 7L9
tél.: 1-888-449-3658
Email: info@novatechmedical.com

ENTHERMICS
Medical Systems
www.enthermics.com
An ISO 13485:2003 certified company

W164 N9221 Water Street • PO Box 443
Menomonee Falls, WI 53052-0443
Phone 262.251.8356 • 800.862.9276
Fax 262.251.7067
Email generalinfo@enthermics.com

ivNow-1

- ➔ **COUNTER TOP UNIT** (cord exits back of unit)
 - ☐ 3L Bag Tilt Kit [5014126]
 - ☐ Pedestal Mount & Bracket Kit [5012245]*
- ➔ **UNIT FOR MOUNTING ON AN EQUIPMENT STAND; STAND INCLUDED** (cord exits bottom of unit)
 - ☐ GCX Light Duty Roll Stand for Devices & Bracket Kit (MOUNT TO GCX EQUIPMENT STAND; STAND INCLUDED) [5012246]
- ➔ **UNIT FOR MOUNTING ON AN EQUIPMENT STAND; STAND NOT INCLUDED** (cord exits bottom of unit)
 - ☐ Equipment Stand Bracket Kit (MOUNT TO EQUIPMENT STAND OR MAST; STAND NOT INCLUDED) [5012735]
- ➔ **WALL MOUNTED UNIT** (cord exits bottom of unit)
 - ☐ 3L Bag Tilt Kit [5014129]
 - ☐ Horizontal Wall Mount Bracket Kit (ONE [1] GCX HORIZONTAL RAIL INCLUDED) [5012734]
 - ☐ Vertical Wall Channel Mount & Bracket Kit (ONE [1] GCX VERTICAL RAIL INCLUDED) [5012241]*
 - ☐ Vertical Wall Channel Mount Bracket Kit (GCX VERTICAL RAIL NOT INCLUDED) [5012288]*

ivNow-2

- ➔ **COUNTER TOP UNIT** (cord exits back of unit)
 - ☐ 3L Bag Tilt Kit [5014127]
 - ☐ Pedestal Mount & Bracket Kit [5012245]*
- ➔ **WALL MOUNTED UNIT** (cord exits bottom of unit)
 - ☐ 3L Bag Tilt Kit [5014130]
 - ☐ Horizontal Wall Mount Bracket Kit (ONE [1] GCX HORIZONTAL RAIL INCLUDED) [5012734]
 - ☐ Vertical Wall Channel Mount & Bracket Kit (ONE [1] GCX VERTICAL RAIL INCLUDED) [5012241]*
 - ☐ Vertical Wall Channel Mount Bracket Kit (GCX VERTICAL RAIL NOT INCLUDED) [5012288]*
 - ☐ GCX Heavy Duty Roll Stand for Devices & Bracket Kit (MOUNT TO GCX EQUIPMENT STAND; STAND INCLUDED) [5013312]


ivNow-3

- ➔ **COUNTER TOP UNIT** (cord exits back of unit)
 - ☐ 3L Bag Tilt Kit [5014128]
 - ☐ Pedestal Mount & Bracket Kit [5012245]*
 - ☐ Pedestal Mount Bracket Kit (PEDESTAL NOT INCLUDED) [5012287]*
- ➔ **WALL MOUNTED UNIT** (cord exits bottom of unit)
 - ☐ 3L Bag Tilt Kit [5014131]
 - ☐ Horizontal Wall Mount Bracket Kit (TWO [2] GCX HORIZONTAL RAILS INCLUDED) [5012242]*
 - ☐ Vertical Wall Channel Mount & Bracket Kit (ONE [1] GCX VERTICAL RAIL INCLUDED) [5012241]*
 - ☐ Vertical Wall Channel Mount Bracket Kit (GCX VERTICAL RAIL NOT INCLUDED) [5012288]*
 - ☐ GCX Heavy Duty Roll Stand for Devices & Bracket Kit (MOUNT TO GCX EQUIPMENT STAND; STAND INCLUDED) [5013312]


ivNow-4 THROUGH -6

- ➔ **WALL MOUNTED UNIT** (cord exits bottom of unit)
 - ☐ Horizontal Wall Mount Bracket Kit (TWO [2] GCX HORIZONTAL RAILS INCLUDED) [5012242]*
 - ☐ Vertical Wall Channel Mount & Bracket Kit (ONE [1] GCX VERTICAL RAIL INCLUDED) [5012241]*
 - ☐ Vertical Wall Channel Mount Bracket Kit (GCX VERTICAL RAIL NOT INCLUDED) [5012288]*


* See equipment diagram below

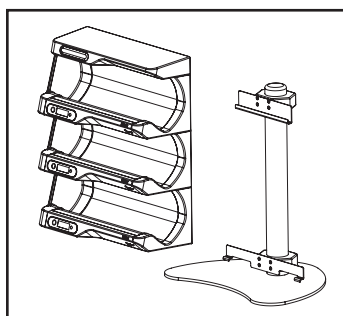


Pedestal Mount & Bracket Kit
5012245
(Brackets only 5012287)

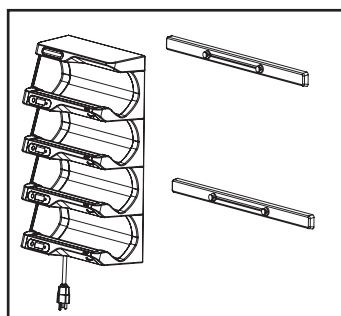

Horizontal Wall Mount Bracket
Kit 5012242 (ivNow-3 to -6)

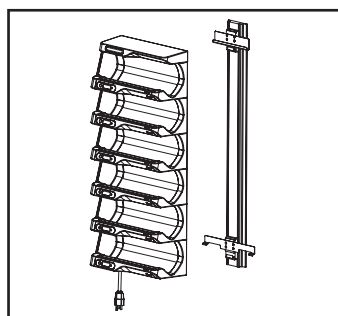

Vertical Wall Channel Mount Bracket Kit 5012241
Brackets only 5012288
(ivNow-6 shown above)

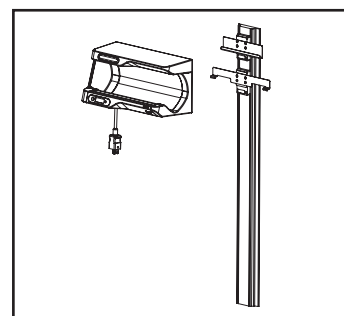

(ivNow-1 shown above)

ENTHERMICS
Medical Systems
www.enthermics.com
An ISO 13485:2003 certified company

W164 N9221 Water Street • PO Box 443
Menomonee Falls, WI 53052-0443
Phone 262.251.8356 • 800.862.9276
Fax 262.251.7067
Email generalinfo@enthermics.com

NOVATECH
MEDICAL

1380 Joliot-Curie, bureau 808
Boucherville (Québec) J4B 7L9
tél.: 1-888-449-3658
Email: info@novatechmedical.com